

MAGELLAN • MAGELLAN BASIC

(ordering code 01.501.042 • 01.501.051) Multifunctional hydrotherapy tub

MAGELLAN - Multifunctional hydrotherapy tub is used to automatic underwater massage, pearl bath and chromotherapy.

MAGELLAN can be used to automatic underwater massage. Electronic control panel and special configuration of 34 jets (5 sections: feet, ankle, knees, hips and spine) let for individual programming working sections. Additionally user can release section which he doesn't want to use during treatment. MAGELLAN multifunctional hydrotherapy tub is equipped with electronic panel lets for programming total time of treatment, time of working particular section, automatic switch off after treatment and enables saving 20 individual programs.

Optionally **MAGELLAN** can be equipped with coloured touch control panel, which allows for programming:

- Total time of treatment
- Time of working every section
- Automatic turn off after finish of treatment
- Free programs - user can save 100 programs
- AutoTest - automatic test of safety working fuses



CHROMOTHERAPY SYSTEM

Emission 7 basic colors and mixed colors

SPOLKA Z O.O.
TECHNOMEX

MAIN FEATURES

The biggest advantage of MAGELLAN tub:

- economic use of water
- not using staff during whole time of treatment
- adjustable foot support for patient
- modern desing and ergonomic construction
- spoutless water filling system
- AutoTest - automatic test of safety working of fuses
- electronic, coloured touch panel - model MAGELLAN

Multifunctional bath tub included the newest system of:

- chemical disinfection - electronic controlled
- automatic release of water - electronic controlled
- release water from installation and pumps - electronic controlled
- automatic filling of tub
- coloured, touch control panel

We offer MAGELLAN in two versions:

- **MAGELLAN BASIC** (electronic control panel)
- **MAGELLAN** (coloured, touch control panel)



Electronic control panel - MAGELLAN BASIC



Coloured, touch control panel - MAGELLAN




STANDARD AVAILABLE COLOURS

	5005 RAL Signal Blue
	5012 RAL Light Blue
	6027 RAL Light Green
	6018 RAL Yellow Green
	6002 RAL Leaf Green
	10217 RAL White
	9001 RAL Cream
	3017 RAL Rose
	2004 RAL Pure Orange

Specify the colour inside and outside:

- Standard colour is blue (5012 RAL) inside.
- Standard colour is white (10217 RAL) outside.

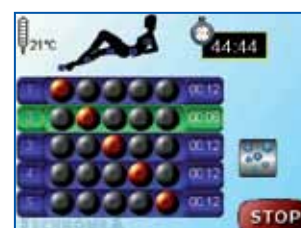

Simple and readable software - MAGELLAN

CE 0197



Electronic, control system allows for connecting:

- underwater massage with pearl bath system (14x5 jets).

MAGELLAN bath tub can be equipped with chromotherapy system - 24 multicoloured diodes (7 main colours and derivated colours).

Independent Digital panel lets for:

- control of time of lighting of colours, individual programming of cycles of alternative lighting.



[1]

SPECIFICATIONS			ELEMENTS OF EQUIPMENT		
Type:	MAGELLAN	MAGELLAN BASIC	Elements of equipment	MAGELLAN	MAGELLAN BASIC
Model:	01.501.042	01.501.051	Hand shower	+	+
Working capacity (I) Level I:	110	110	Automatic underwater massage - 34 jets - 5 sections (foots, ankle, knees, hips and spine)	+	+
Working capacity (I) Level II:	220	220	Pearl bath system - 14x5 jets	+	+
External Dimensions (mm)(LxWxH):	2200 x 950 x 980*		Filling system	automatic	automatic
			Setting of desired water level	automatic	automatic
			Automatic release of water - Electronic controlled	+	+
*Optionally the tub is available with the height of 70 cm. STANDARD EQUIPMENT • Stair 53L x 43W x 39H cm [1] - (19.037.001)			Control panel	coloured, touch	electronic
			AutoTest - automatic test of safety working fuses	+	-
			Pump "dry-working" protection	+	+
			Timer	+	+
			Automatic switch off after time of treatment	+	+
			Chromotherapy	+	option
			Disinfection system	+	option
			Disinfection system with ultraviolet light	option	option
			Electrotherapeutic bath system	option	option
			Aromatherapy	option	option
			System of keeping temperature/heating up the water in the tub	option	option