





Patient database

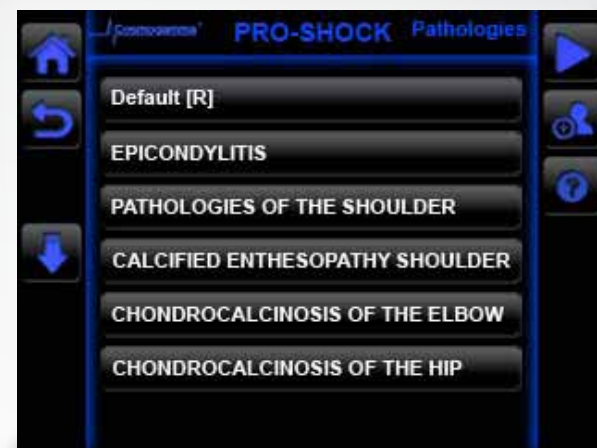

Pathologies

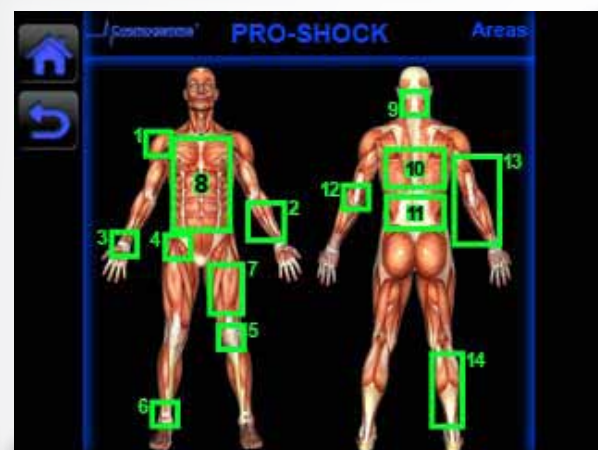

Anatomical areas



Emission treatment



PRO-SHOCK Waves

THE MOST VERSATILE SHOCKWAVES EVER

- **EXTREMELY REDUCED ENERGY** TRANSMITTED TO OPERATOR
- **EXTREMELY SIMPLE** TO OPERATE
- HUGE QUANTITY OF **PRESET PROGRAMS**
- POSSIBILITY TO **CREATE/SAVING** YOUR OWN PROGRAMS
- **E.S.A. SYSTEM** (EXTREMELY SOFT ABSORBED SYSTEM)



AC International s.r.l.

via Canapa 22 - 44042 Cento (FE) - Italy - Tel. +39 051 6836284 - Fax. +39 051 6831061 - www.cosmogamma.com - info@acintl.it

PRO-SHOCK WAVES

(CYB.SHOCK) SHOCK WAVES THERAPY

Pro-Shock Waves incorporates a computer with "Windows CE" system that provides almost infinite programming capabilities and allows the operator to enter the patient's data including information on the history and evolution of patient's disease.

A convenient external **USB** allows transfer of patient's data to an external database, as well as software updating.

The use of modern mechanical technologies applied to this machine guarantee the best quality of delivery of the shock waves on the patient. Furthermore, there are several energy transmitters with different focalizations, adopting the use of the machine to different types of application on the human body.

From sports injuries to chronic diseases and disorders requiring rehabilitation, **Pro-Shock Waves** offers the ultimate technology in shockwaves and, with a wide range of pre-set protocols, provides highly effective results in the treatment of pain.

In addition to the large amount of pre-treatment programs, through visualization of different anatomical areas, you can select all the treatments dedicated to a specific area of the human body.

The unique patented **E.S.A.® System (Extremely Soft Absorbed)** reduces the energy transmitted to the therapist up to 90%. With its innovative and patented System, the Pro-Shock Waves provides unequalled parameter flexibility, allowing the therapist to reduce the energy and create the perfect combinations of parameters for each patient.

The **Pro-Shock Waves** has a unique, sophisticated design that ensures maximum functionality for a variety of therapeutic applications.

SPECIFICATIONS:

Pro-Shock Waves is a very advanced product delivering radial shockwaves for therapeutic use. It is absolutely one of the most versatile devices on the market:

- Selection by anatomical area of the body
- Selection by medical branches
- Pulsed mode up to 20 Hz
- Pressure power up to 5 bar
- Possibility to create custom programs

- Patient Database
- Huge quantity of preset programs
- E.S.A.® System - This exclusive system reduces more than 90% of energy transmitted to the therapist during the treatment. This system works as a classic shock absorber and removes any problems from the transfer of energy from the handle to the therapist.

TECHNICAL DATA:

Number of output:	1	Handpiece interchangeable with transmitters.
Power supply:	50-60 Hz (±10%)	part interchangeable supplied
- 230 Vac	50-60 Hz (±10%)	
- 115 Vac	50-60 Hz (±10%)	
Maximum power input from the mains:	80 VA	
Double fuse protection on the mains:		
- 230 Vac	2x2 A-T - 5x20 mm	
- 115 Vac	2x4 A-T - 5x20 mm	
Indications:	Backlit graphic LCD display 320 x 340 pixels Touch Screen	
Type of emission of the pulses:	Single, Continuous, Burst	1 pc. - 15 mm for
Type of emission Burst:	(1÷100) burst pulses (10÷2000) ms pause between burst pulses	1 pc. - 15 mm m
		1 pc. - 9 mm mu
		1 pc. - Gel tube
		1 pc. - 0.5 bar
		1 pc. - 1 Hz
		• 36 mm trans
		stage
		rent
Weight:	7.2 kg	
Dimensions (LxWxH):	40 x 42 x 110 cm	
Safety class:	I type BF	
Protection degree:	IP20	

STANDARD

OPTIONAL

